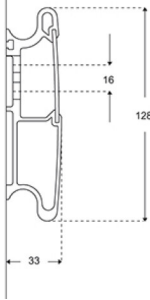




PRESSALIT®

Product	PLUS wall track, 900 mm
Item number	R9809112
Colour	<p>VVS:615841109 NRF:7501014 EAN:5708590305752 HMI: 71740</p> <p>Anthracite grey</p>
Product description	 <p>PLUS wall tracks are used for mounting of support arms, shower seats and wash basin brackets, and allows sideways adjustment of these products. See below a selection of products for fitting to the wall track.</p> <p>PLUS wall track, 900 mm, for horizontal mounting. The cover plate hides bolts and wires and facilitates cleaning. Materials: Aluminium and polystyrene. See below a selection of products for fitting to the wall track. Excl. mounting kit.</p>
Adjustability	Fixed mounting
Product capacity	Maximum load: Up to 1000 kg depending on the wall type.
Materials	The wall track is primarily made of 10-micrometre anodized aluminium. The cover plate is made of polystyrene and the ends of the track are capped with ABS polymer.
Ergonomy	The track distributes the weight of the individual bathroom product, for example a wash basin bracket, along its entire length. The products can be adjusted and locked on to the track with intuitive locking functions.
Features	The horizontal track makes it possible to shift the product along the wall to constantly optimize the available space in the bathroom. The front cover plate conceals mounting bolts and any cables.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



PRESSALIT®

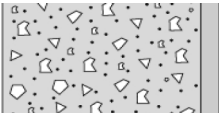
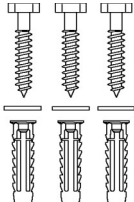

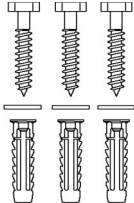
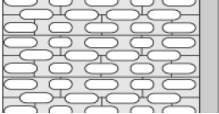
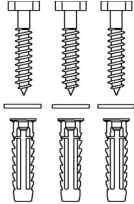
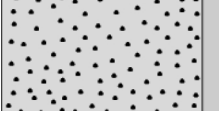
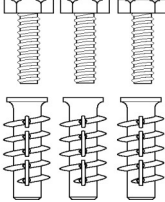
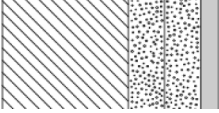
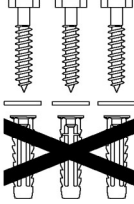
Product	PLUS wall track, 900 mm
Item number	R9809112
Cleaning	<p>PLUS wall tracks.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
Declaration of conformity	<p>MDR-ACC This product fulfils all relevant requirements in Regulation (EU) 2017/745 on medical devices in risk class I. This product is considered to be an accessory for a medical device, as it is not itself a medical device, but is intended to be used together with several particular medical devices to specifically enable the medical devices to be used in accordance with their intended purposes.</p>



PRESSALIT®

Product	PLUS wall track, 900 mm
Item number	R9809112

Installation in different wall types

	<p>Concrete >B25</p> <p>V8620 - Mounting kit, 3 coach screws Ø8 x 70, plugs SX10. To be ordered separately.</p>	
	<p>Solid Brick >= Mz12</p> <p>V8620 - Mounting kit, 3 coach screws Ø8 x 70, plugs SX10. To be ordered separately.</p>	
	<p>Solid Sand Lime Brick >= Hlz12</p> <p>V8620 - Mounting kit, 3 coach screws Ø8 x 70, plugs SX10. To be ordered separately.</p>	
	<p>Aerated Concrete P4</p> <p>V8670 - Mounting kit, 3 bolts M8 x 35, plugs FTP M 8. To be ordered separately.</p>	
	<p>Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood</p> <p>V8620 - Mounting kit, 3 coach screws Ø8 x 70, plugs SX10. NB! Plugs are not needed. To be ordered separately.</p>	



PRESSALIT®

Product	PLUS wall track, 900 mm
Item number	R9809112

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®